

The next
generation in
capsule sorting
from MOCON


VERICAP[®] 4000

Highest

Throughput

with High

Sorting Precision

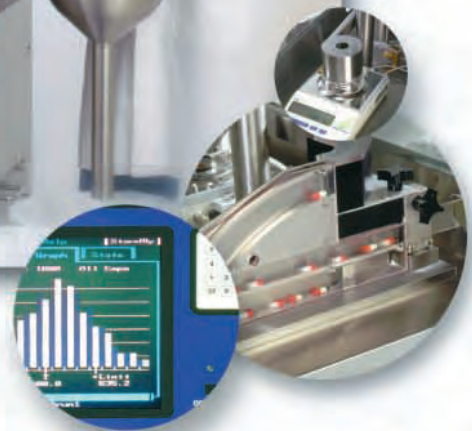

High-Speed Capsule
Weighing and Sorting System



Equipos de control de calidad

Tel.: (+34) 902 450 160

Fax: (+34) 902 433 088

ermec@ermec.net


www.ermec.net

mocon[®]

The VERICAP[®] 4000 features...



The Fastest Capsule Weight Sorter in the World

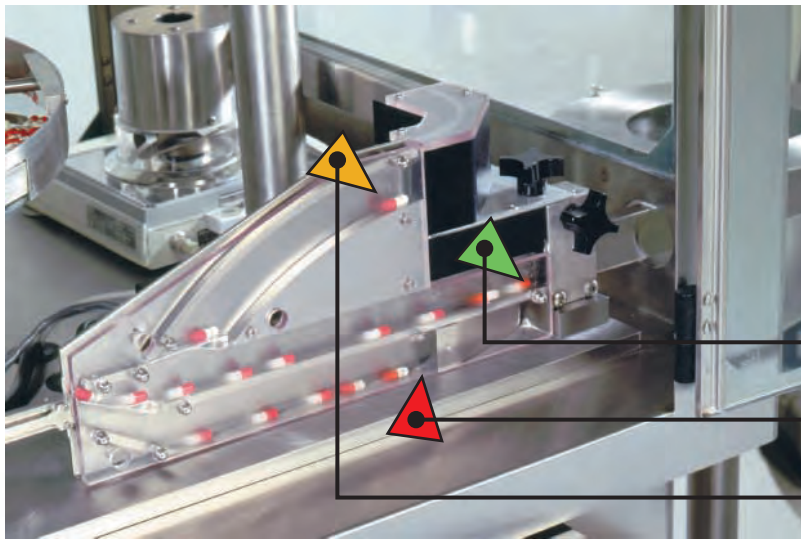

The high-speed sensor featured in the VERICAP[®] 4000 is capable of weighing capsules at speeds as high as 132,000 capsules per hour. To ensure precision, the VERICAP 4000 sensor re-zeros itself between every capsule weight measurement.

- **“Accepted” Capsules**
- **“Rejected” Capsules**
- **“Checkweighed” Capsules**


System Enhancements:

- Improved two flipper system
- New data archival with Zip[®] drive
- Automatic gain adjustments
- Low-electronic-noise power supplies and motors
- Greater software precision
- Increased use of stainless steel throughout
- Quick disconnects for easier cleaning
- New glass doors
- Checkweigher included standard
- Printer included standard

Increased Accuracy

The VERICAP 4000 offers increased accuracy over previous models.

This increased accuracy is the result of a number of proprietary enhancements including improved electronics and greater software precision.

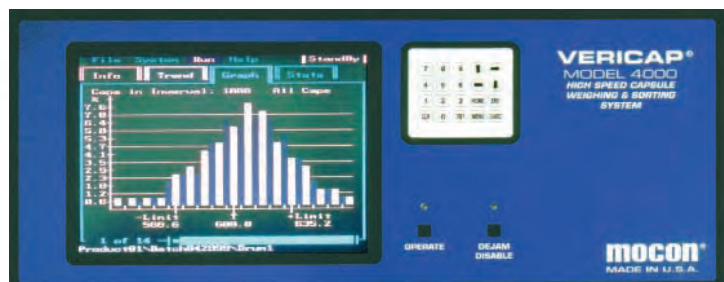

Throughput

The VERICAP 4000 is easy to use, requiring the calibration of **only one sensor!** Initial set-up and calibration is much easier than competitive systems. It also features an integrated, on-board checkweigher, which allows the VERICAP 4000 to calibrate itself **automatically** during operation.

Up to 132,000 capsules per hour



Integrated On-Board Checkweigher



The VERICAP 4000 checkweigher provides a built-in analytical balance for a gravimetric weight check of a capsule. This capability automatically makes minor adjustments to the VERICAP 4000 system calibration should it become necessary during the run. This ensures maximum precision and traceability. These checkweighed capsules are collected separately and are not returned to the batch population. The frequency and number of capsules checked is adjustable.

Simple and Supported Validation



On-site assistance by MOCON personnel

The VERICAP 4000 features several available validation programs based on the proven VERICAP 2500 validation packages. Depending on customer requirements, the VERICAP 4000 comes with optional validation assistance package or full validation.



Equipos de control de calidad

Tel.: (+34) 902 450 160

Fax: (+34) 902 433 088


ermec@ermec.net

www.ermec.net

Easier to Maintain

The VERICAP 4000 features increased use of stainless steel and glass for durability and cleanability. In addition the system has quick disconnects for all table top air lines and electrical connections to further enhance cleanability and reduce machine tear down and re-assembly time.

Capsule Flow Diagram

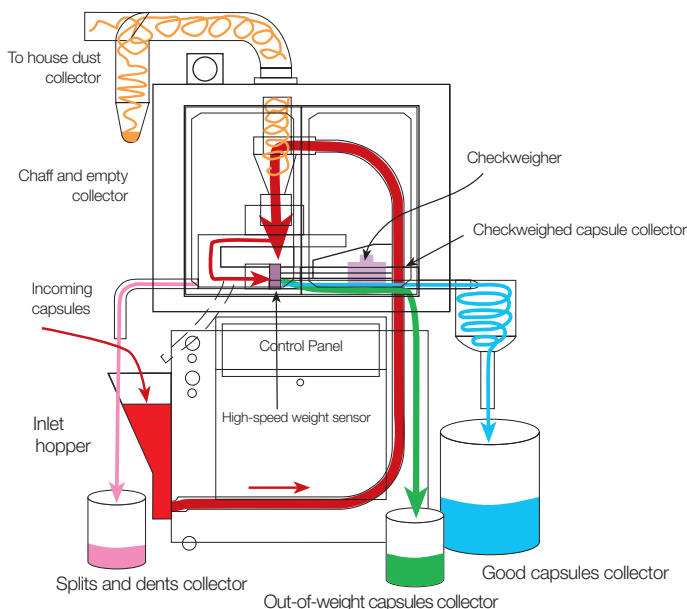

- Incoming capsules from encapsulator or drums. █
- Empties, chaff, halves, are removed and collected. █
- Splits, dents, and other jams are cleared and collected. █
- Out-of-weight tolerance capsules are rejected and collected. █
- "Good" capsules are collected. █

Capsule Change Parts



Additional change parts are available for size #00, 0, 1, 2, 3 and 4 capsules. Changeover is easy and fast, without the need for tools. One change part is included in the purchase price of the VERICAP.

Fast Payback in Owner Surveys

Fast payback times have been experienced by many VERICAP users. Some companies have paid for the machine in just one batch of product. If desired, MOCON would be happy to prepare a payback analysis based on your specific capsule operation.

Specifications

- Speed: Up to 132,000/hr (2,200/minute)
- Inlet Hopper: Volume is 1033 in³ (16.9 liters)
- Power: 120 VAC±10% @ 60 Hz, 1.5A
All international voltages available
- Noise Level: Less than 75 dbA with cover closed.
- Air Requirements: Main unit 100 psi (7 bars).
Inlet supplied with quick disconnects to 3/8" (0.95 cm) hose.
- Vacuum Requirements: Soft vacuum for dust collection, 4" (10.2 cm) duct.

Note: For maximum throughput, room humidity should be maintained between 40 and 60% RH.



Tel.: (+34) 902 450 160
Fax: (+34) 902 433 088
ermec@ermec.net
www.ermec.net

mocon[®] Phone 612-493-6370
Fax 612-493-6358
7500 Boone Ave N, Minneapolis, MN 55428

Visit our website at
www.mocon.com
E-mail: mocon@mocon.com

MOCON and VERICAP are registered trademarks of MOCON, Inc.
Zip is a registered trademark of Iomega Corporation
Copyright © 1999 MOCON, Inc.
MOCON reserves the right to change specifications without notice as part of our continuous program of product improvement.